

WARRANTY

To activate the **lifetime** warranty, please send the order ID to our official designated email address:

warranty@anzeser.com

Contact us

support@anzeser.com

Contact us through above official email address to get faster customer support.

EC

REP

eVatmaster Consulting GmbH
Bettinastr.30
60325 Frankfurt am Main, Germany
contact@evatmaster.com

UK

AR

UK CROSSBORDER LIMITED
7 Bell Yard London WC2A 2JR, UK
United Kingdom

Magnetic Stirrer with Heater



Lab Supplies

This hot plate contains an integrated magnetic stirrer, ideal for use in a laboratory. The dual controls allow you to set heat output and stirring speed independently. The hot plate can reach a maximum temperature of 350 °C (662 °F) with an output of 180 watts. The magnetic stirrer is adjustable from 100 to 1600 rpm.

A support stand with an adjustable thermometer holder is also included. The support stand extends about 10" above the hot plate's surface.

Components:

Base with power switch, two knobs and attached hot plate 1" Teflon coated magnetic stir bar

Power cord

Support stand (base rod, extender rod, right angle clamp and thermometer holder)

You may also want the following (not included):

Partial immersion or digital thermometer

Glass beaker or flask

Operation

To set up the support stand: Screw the two rods together making one support rod. Make sure they fit tightly. Insert the slotted end of the rod into the plastic holder on the back of the base. Slip the right angle clamp onto the support rod and tighten the knob. Slide in the thermometer holder and tighten the other knob. The height of the thermometer holder can be adjusted by moving the clamp up or down on the support.

To use the hot plate:

1. Make sure the knobs are turned off. The white mark on each knob should be pointing down.
2. Connect the power cord to the back of the base. The base should be set on a hard surface, like a countertop or table.
3. Plug the other end of the power cord into a standard outlet. Make sure the outlet is safe and the power cord is not twisted.
4. Flip the power switch on the right side of the base. The line (-) indicates on and the circle (o) indicates off.

WARNING: Be very careful with the hotplate turned on. Do NOT touch the aluminium surface as it gets very hot and can cause severe burns.

5. Rotate the right knob in a clockwise direction. At maximum power the hot

plate will take around 20-30 minutes to boil a large beaker of water. The surface in the middle of the plate can reach 350°C (662°F) when the knob is turned all the way to maximum power.

6. Adjust the knob to keep the desired heat level, using a thermometer to check the temperature of the solution.
7. To turn off the hot plate, rotate the knob counter-clockwise until you hear a click and the white mark is pointing down.
8. Make sure both knobs are off and flip the power switch off. Unplug the cord from the outlet by pulling on the plug, not the cord.

To use the magnetic stirrer:

1. Drop the stir bar into the solution and let it settle.
2. Turn on the magnetic stirrer by flipping the power switch on (-), then rotating the left knob in a clockwise direction.
3. The further you rotate the knob, the faster it will spin to stir your solution.
4. To turn off the magnetic stirrer, rotate the knob counter-clockwise until you hear a click and the white mark is pointing down. Flip the power switch off. Unplug the device.
5. Remove the stir bar with a metal utensil. The stir bar will stick to it so you can lift it out easily.

Main Technical Parameters:

Model	Voltage (V)	Speed (r.p.m)	Heating Power(w)	Max. Volume (ml)	Max.Temp (plate)	Plate size (mm)
SH-1	☒ 110V	100~1600	—	2000	350℃ ± 10%	120×120
SH-2	☐ 220V		180			3000
SH-3	☐ 230V		500	5000		
SH-4	☐ 240V		600			

**Final Details for Order #114-6482155-3692218**[Print this page for your records.](#)**Order Placed:** September 17, 2023**Amazon.com order number:** 114-6482155-3692218**Order Total:** \$59.54**Shipped on September 18, 2023****Items Ordered****Price**

1 of: ANZESER Magnetic Stirrer Hot Plate w/Stir Bar, 100-2000rpm Magnetic Hotplate Stirrer, 200W Heating Power Hot Plate Stirrer, Max 520°F Hot Plate with Magnetic Stirrer \$55.00

Sold by: Best Magnetic Stirrer Store ([seller profile](#)) | Product question? [Ask Seller](#)

Supplied by: Best Magnetic Stirrer Store ([seller profile](#))

Condition: New

Shipping Address:

David Malawey
3500 GREEN RIDGE CIR
BRYAN, TX 77802-4715
United States

Shipping Speed:

Delivery in fewer trips to your address

Payment information**Payment Method:**

Visa ending in 4670

Item(s) Subtotal: \$55.00

Shipping & Handling: \$0.00

Billing address

David Malawey
3500 GREEN RIDGE CIR
BRYAN, TX 77802-4715
United States

Total before tax: \$55.00

Estimated tax to be collected: \$4.54


Grand Total: \$59.54

Credit Card transactions

Visa ending in 4670: September 18, 2023: \$59.54

To view the status of your order, return to [Order Summary](#).

[Conditions of Use](#) | [Privacy Notice](#) © 1996-2023, Amazon.com, Inc. or its affiliates

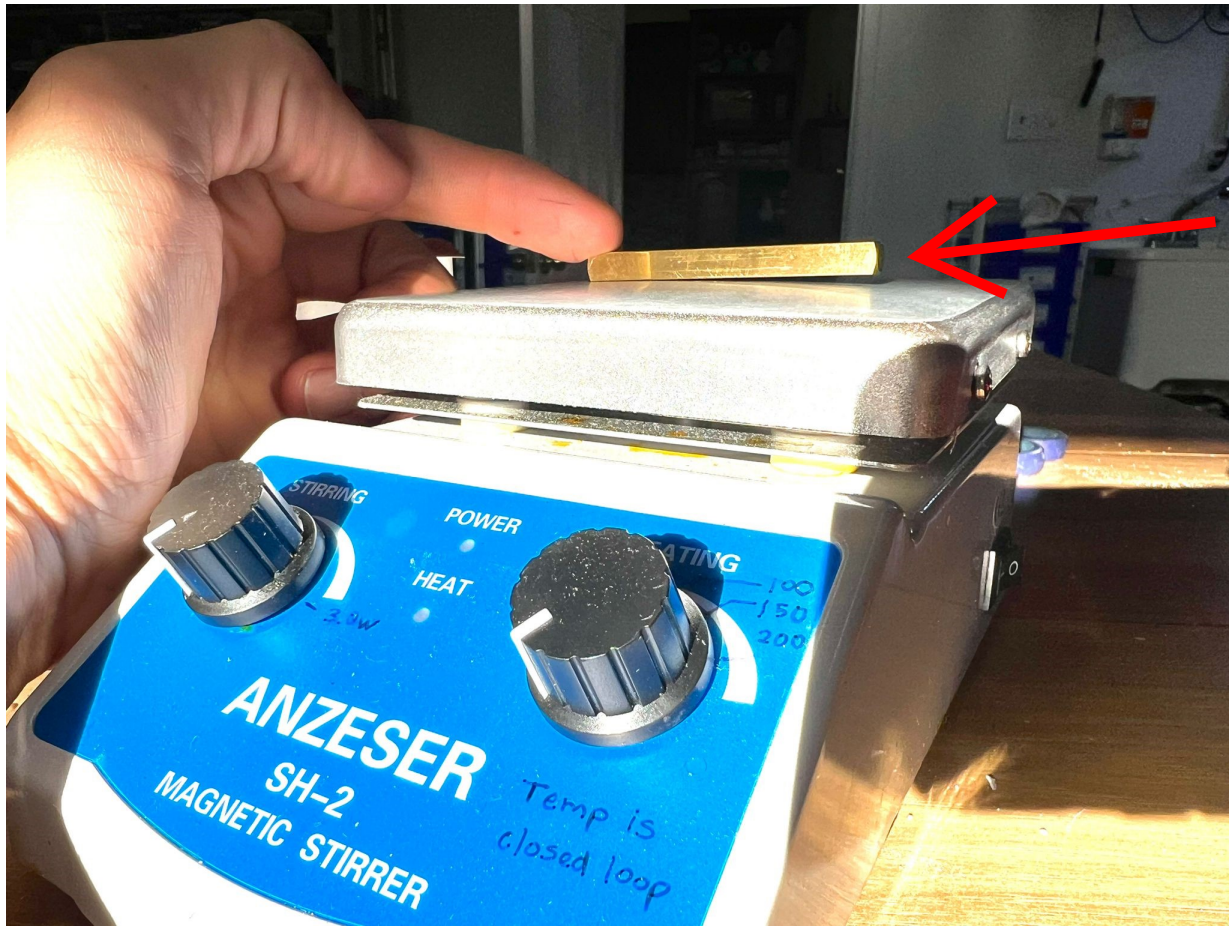

Issue: The top of the plate is not flat. When the product was purchased, it was slightly warped and now after some heating cycles it is warped strongly. glasses cannot sit flat on the top and have a steady operation. - David, 2023.12.12